

RASCAL[®]

VENTURE BY MEGA MOTION[®]



AVAILABLE IN BOTH 3
AND 4-WHEEL MODELS

Width adjustable
armrests

Rear storage
pocket

17" (W) x 17" (D)
compact, foldable,
height-adjustable
seat with 90° swivel

Swift Release Technology
for quick disassembly

Lightweight 12 Ah
SLA battery pack

Easily accessible
charging port

Front basket
included

Easy-to-grip tiller
adjustment knob

Non-scuffing,
flat-free tires

Color-through
shroud

FDA Class II Medical Device *




LONG LASTING
RANGE UP TO
10 MILES
AT 200 LBS.
ON 3-WHEEL MODEL



4 MPH
TOP SPEED



**TIGHT
TURNS**

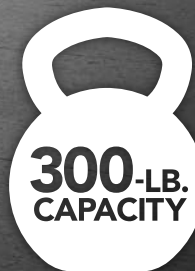

**300-LB.
CAPACITY**



RASCAL[®]
POWER MOBILITY

rascalpowermobility.com

* Rascal[®] Venture Scooters are FDA Class II Medical Devices designed to aid individuals with mobility impairments.



SWIFT RELEASE TECHNOLOGY provides a quick way to transport your Venture scooter in the trunk of your vehicle in just a few simple steps.



Rascal® Venture is available in both 3 and 4-wheel models.



Intuitive controls include horn, speed adjustment, battery meter, key switch, and forward/reverse throttle.



Recharging the Venture is fast and easy with the charging port located conveniently below the right handlebar.



Maximize your comfort with flip-up armrests, adjustable tiller, and seat height adjustment.

SPECIFICATIONS

WEIGHT CAPACITY	300 lbs.
GROUND CLEARANCE	3-wheel: 1.26" 4-wheel: 1.36"
MAXIMUM SPEED¹	4 mph
TURNING RADIUS²	3-wheel: 33" 4-wheel: 46.75"
LENGTH²	3-wheel: 38.5" 4-wheel: 40.5"
WIDTH²	22.5"
HEIGHT	34.75"
TIRES	8" solid (flat-free)
RANGE PER CHARGE^{1,3}	3-wheel: Up to 10 miles at a user weight of up to 200 lbs. (8 miles at over 200 lbs.) 4-wheel: Up to 8 miles at a user weight of up to 200 lbs. (6.6 miles at over 200 lbs.)
TOTAL WEIGHT	3-wheel: 96.5 lbs. 4-wheel: 101.5 lbs. (includes battery pack)
SEATING²	Compact, foldable, black vinyl seat with 90° swivel and height adjustment (22.7 lbs.)
DRIVE SYSTEM	Rear-wheel drive, single motor, 24-volt DC motor
DUAL BRAKING SYSTEM	Regenerative and electromechanical
BATTERY REQUIREMENTS	12 Ah sealed lead acid, non-spillable deep-cycle removable battery pack
BATTERY CHARGER	2 amp off-board charger
MAXIMUM RATED SLOPE	6°
MAXIMUM OBSTACLE CLIMBING²	2"

¹ Varies with base model, user weight, terrain type, battery amp hour rating (Ah), battery charge, battery condition, motors, controller type, tire type, and tire condition. This specification is subject to a variance of +10%, -5%.

² Due to manufacturing tolerances and continued product improvement, this specification can be subject to a variance of (+ or -) 3%.

³ Tested in accordance with ANSI/RESNA, WC Vol. 2, Section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

NOTE: This product conforms to all applicable ANSI/RESNA and ISO 7176 Series standards. All specifications subject to change without notice.



RASCAL
POWER MOBILITY

rascalpowermobility.com

* Rascal® Venture Scooters are FDA Class II Medical Devices designed to aid individuals with mobility impairments.